



DR-2000

GENERAL OPERATING TABLE

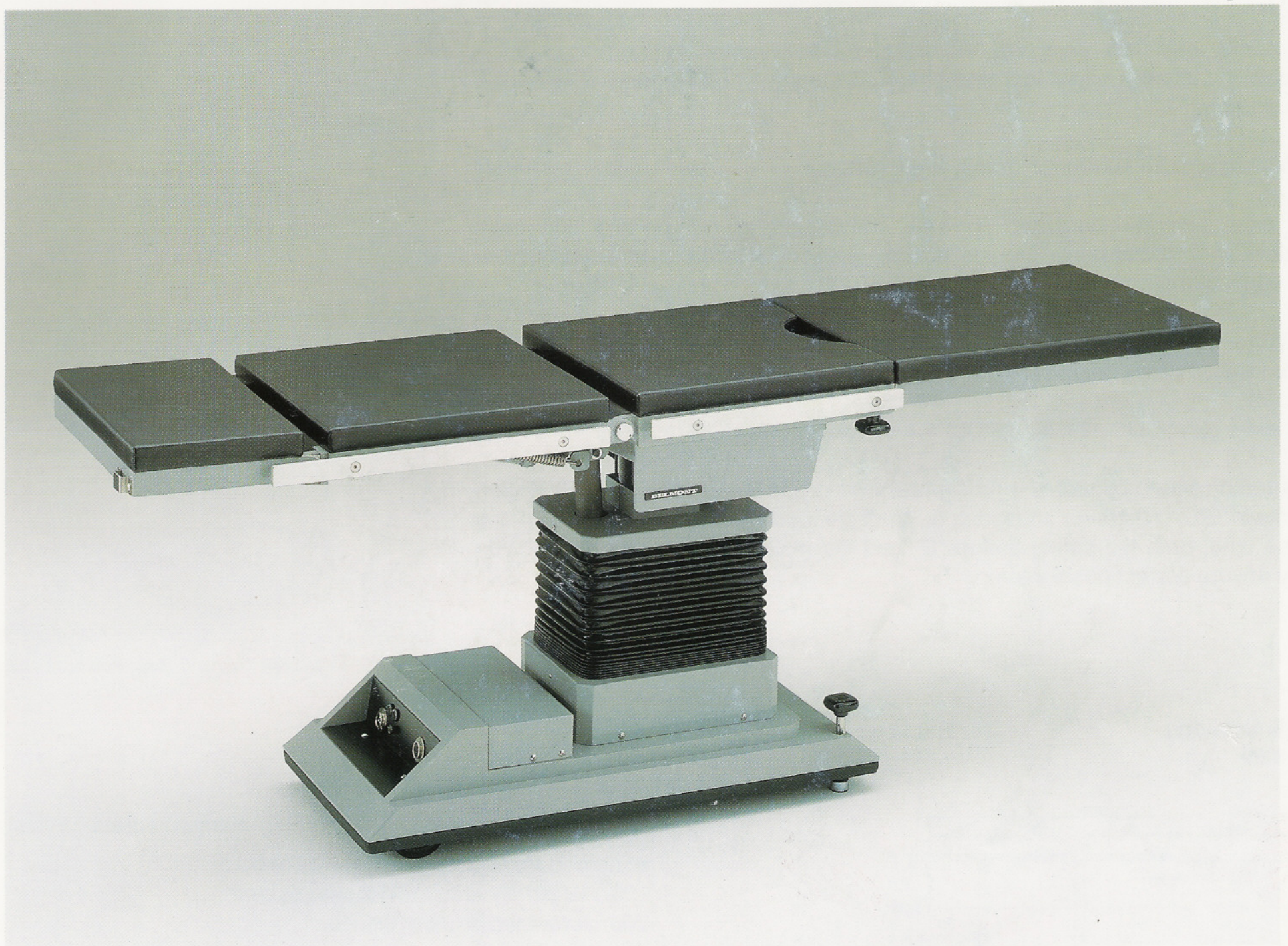

 Belmont®

GENERAL OPERATING TABLE DR-2000

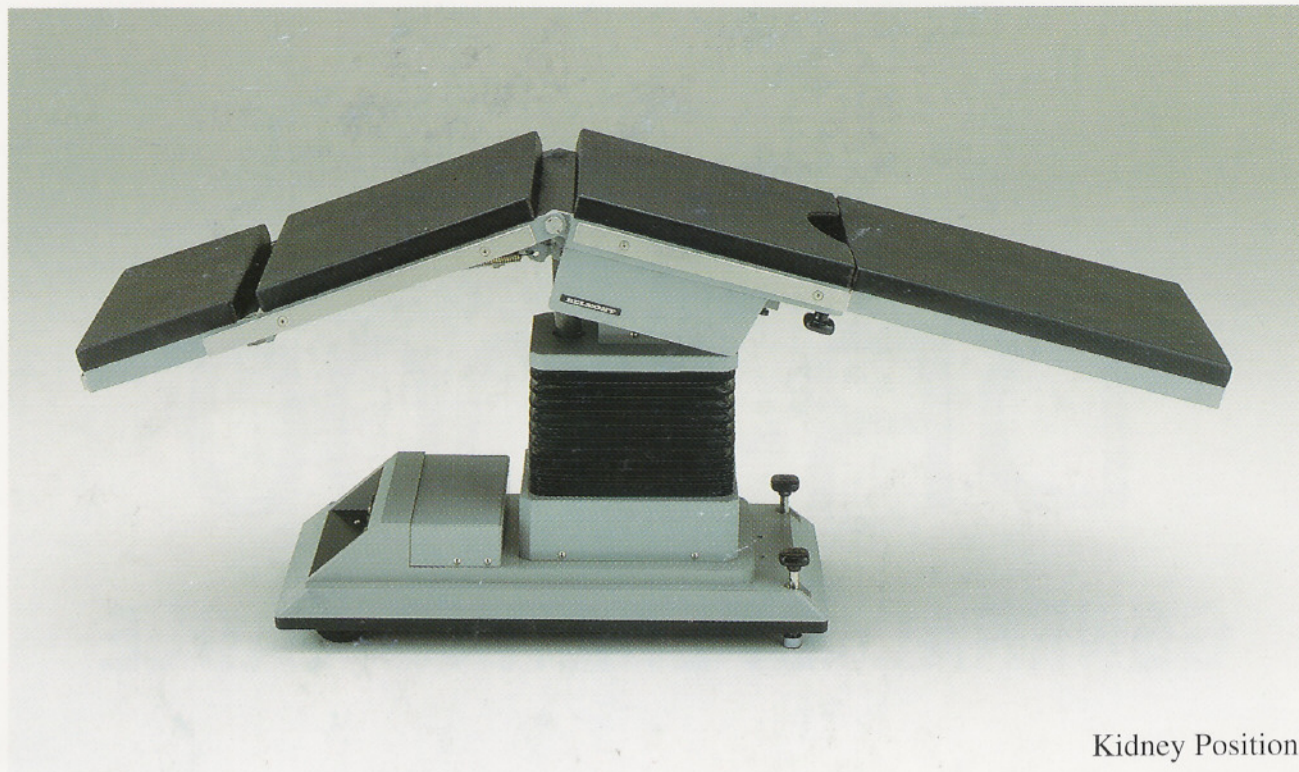

Kidney Position

The DR-2000 is a general purpose operating table providing stability and long life at economical price.

Belmont's own low pressure electro-hydraulic system operates the table smoothly without a slight jerk and shock at the start and the end of hydraulic movements.


The table top lifts to 900mm from a low of 600mm.

20° each Trendelenburg and Reverse Trendelenburg, 50° up and 35° down back support movements provide patient positioning to the best advantage of the surgical team.

An Automatic Leveling Device will bring the table top to a horizontal position from trendelenburg and reverse trendelenburg tilt angles.

For extra safety, DC.12V is used for the control circuit.

The hand remote control is a flush-proof membrane type with the safety button marked "E". Whatever movements shall be obtained by depressing at the same time, both the safety button "E" and the pictographic button so that you may ensure correct positioning.

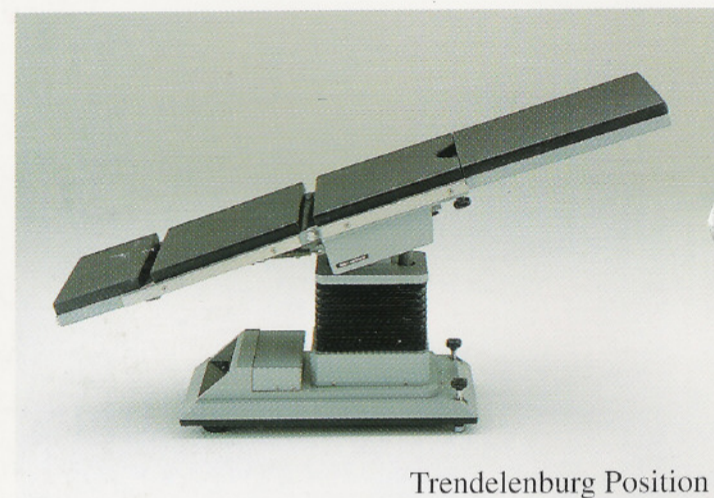



Trendelenburg Position



Gynecological Position

(Douche Tray and Thigh Rests are optional.)

Dimensions -mm-

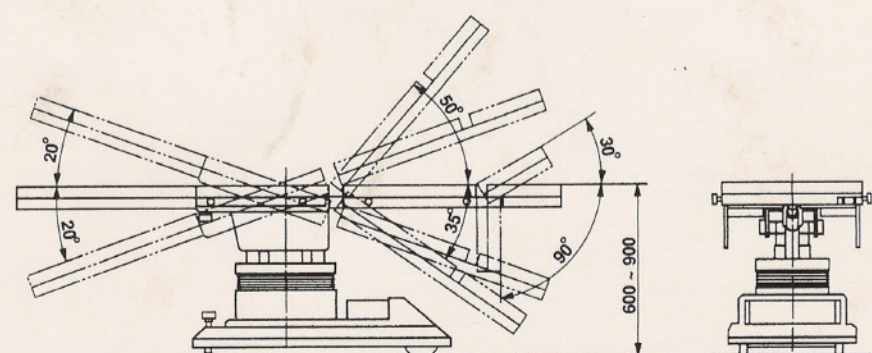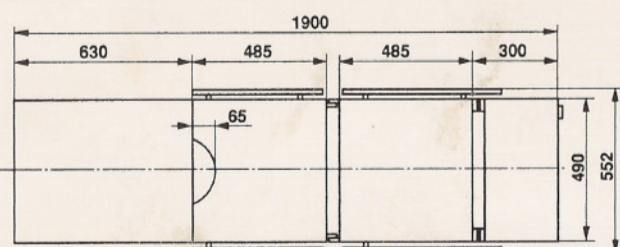

Specifications

Table top	1900×490mm
Table top height	600~900mm By electro-hydraulic system
Trendelenburg and Reverse trendelenburg	-20° ~ +20° By electro-hydraulic system
With automatic leveling device	
Back support angle	-35° ~ +50° By electro-hydraulic system
Headrest angle	90° ~ +30° By manual control
Power voltage	AC. 100/115/220/240V
Control voltage	DC.12V
Net weight	145kg

The manufactures reserve the right to alter and improve the specification without previous notice.

Belmont®

Printed in Japan DR-2000 9503AK